

# Hygiene Monitoring Kit FUKITORIMASTER

## Rapid Protein Detection Kit for Cleanliness Test

**Protein Residue**


**10 minutes**

**Restaurant !**

**4 Level Colors**




### ■ Procedure

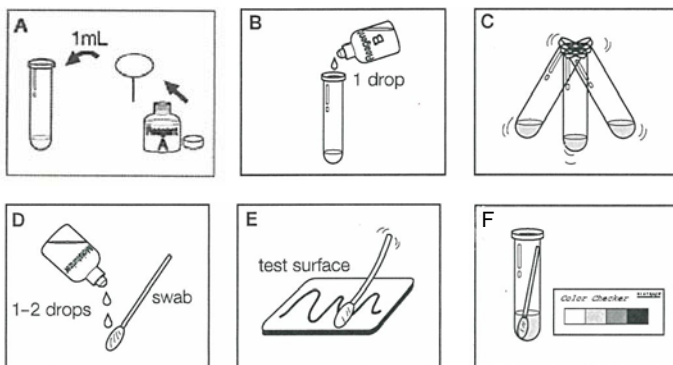

- Pour 1mL of Reagent A by dropper
- And pour 1 or 2 drops of Reagent B
- Mix it well
- Wet a swab with Moisturizer
- Swab on the surface (10cmx10cm)
- Put the swab into tube and wait for 10 minutes and check the solution color by color chart

### ■ High Sensitivity Data (Temperature and Time Parameter)

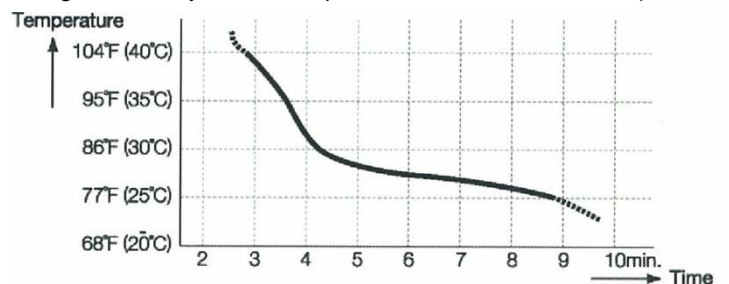

### ■ Contents of Kit

Reagent A (2x100mL), Reagent B (1x20mL), Swab(2x100本)、Tube(5本)、Dropper 1mL(1)、Moisturizer (1x20mL)、Color Checker (1)

**AR BROWN**  
Think Quality & Create Value

Dai-ni Marutaka Bldg.,7-13-8  
Ginza, Chuo-ku, Tokyo 104-0061 Japan

Tel:+81-(0)3-3545-3993 Fax:+81-(0)3-3543-8865  
Email: a-r@arbrown.com

| Code | Product        | Content     | Retail Price |
|------|----------------|-------------|--------------|
| F200 | Fukitorimaster | 200Test/box | ¥27,000      |